


# CXDI-401G COMPACT

UPGRADING TO DIGITAL  
RADIOGRAPHY WITH EASE

*We Speak Image*

The CXDI-401G COMPACT Digital Radiography System with improved resolution allows for easy upgrades to your existing radiography equipment and fits easily into most universal Bucky systems. This unit can be retrofitted into a range of radiography devices, such as upright stands and tables. Streamline your radiographic workflow with the CXDI-401G COMPACT.

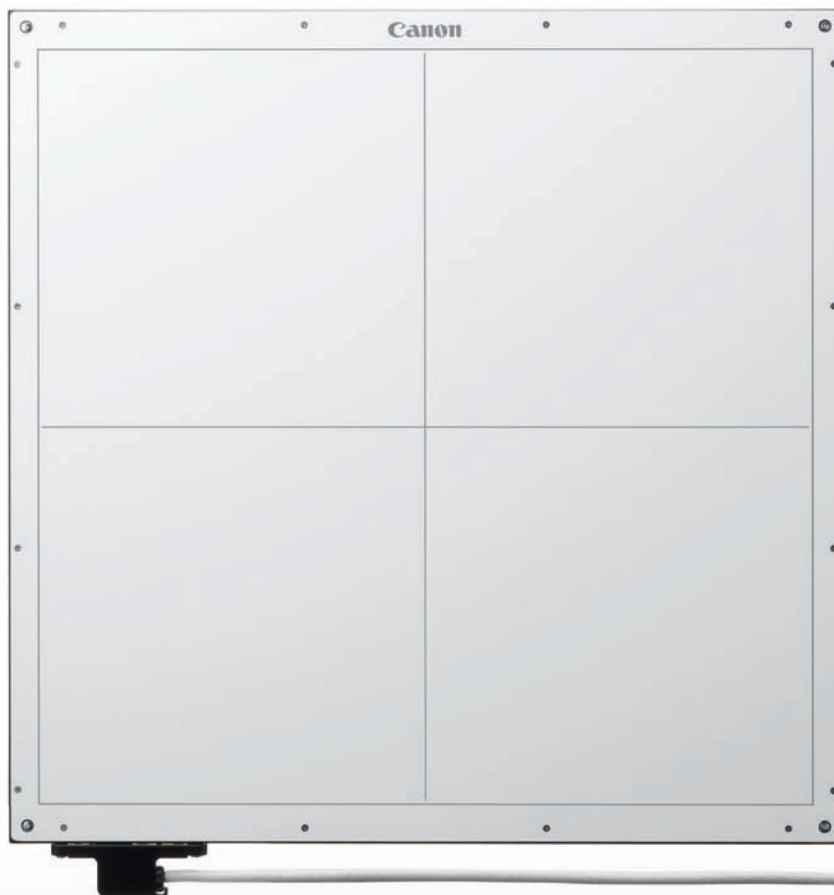

**Canon**

## CXDI-401G COMPACT Specifications

- Fit into 46x46cm Bucky without modification
- Thickness: 15mm
- Grid and AEC are attached to bucky unit

Imaging area	43cm x 42cm
Sensor dimensions	46cm x 46cm x 1.5cm
Pixel pitch	125 µm
Resolution	11.3 mega pixels
Scintillator	GOS
Preview image display	Approx. 3 sec.
Full-processed image display	Approx. 5 sec.
Cycle Time	Approx. 15 sec.
Sleep to ready	Approx. 10 sec.



Supporting CXDI Control Software NE

Up to four detectors can be used in one room, offering you maximum flexibility to accommodate the required workflow in your institute.



**Canon**

Canon Medical Systems Europe  
 Bovenkerkerweg 59-61  
 1185 XB Amstelveen  
 The Netherlands  
[www.canon-europe.com/medical](http://www.canon-europe.com/medical)